



Header


# Header

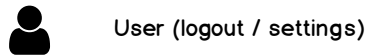


clac : home (Inbox)

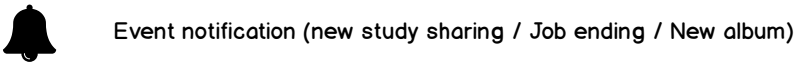

When the user clic on a notification (event or comment)  
He will be redirect to the album, comment or studies

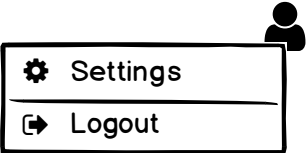


When User clic

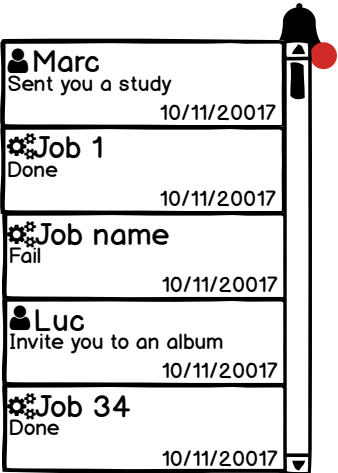

When Event notification clic

Sort by date descending  
Notification source :  
User  
Job  
Notification date  
Notification  
Job  
Done  
Fail  
Album  
Paul invite you to the Album\_name  
Paul add studies to the Album\_name  
Studies  
Paul sent you a studies

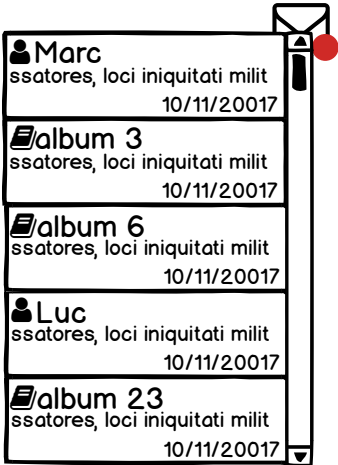

When Message notification clic

Sort by date descending  
Notification source :  
User  
Album  
Notification date  
Notification  
Message



Kheops logo

Tabs bar


 **Inbox**


 Albums


 Favorites


 Recents


 Jobs

 **Inbox** ...

 **Albums** ...

 **Favorites** ...

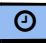 **Recents** ...

 **Jobs** ...

Inbox : show all studies for the connected user


Album : show all albums for the connected user


Favorites : show all the favorites studies for the connected user


Recents : show all recents studies for the connected user


Jobs : show all jobs for the connected user


Tabs bar idea from GitHub :


 **Usirix-Foundation**

 **Repositories** 11

 **People** 6

 **Teams** 0

 **Projects** 1

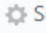 Settings

Type: All ▾

Language: All ▾

Customize pinned repositories

New

**KheopsMockup**  
Updated 10 hours ago

Top languages  
Java Shell JavaScript Lua

# Study list



<input type="checkbox"/>	Patient Name	MRN	Accession #	Study Date	Modality
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 20 2007	CT
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 28 2007	MR
▶	<input type="checkbox"/> ENTERIX	foidh	92485743	Oct 19 2010	PT/SC
▶	<input type="checkbox"/> MACOESSIX	UGVdzg3	987645	Sep 19 1999	PT/SC
▶	<input type="checkbox"/> MANIX	h4873g	0	Nov 1 2003	CT
▶	<input type="checkbox"/> PHENIX	P-000002	A389274983	Mar 21 1996	CT
▶	<input type="checkbox"/> PHENIX	P-000002	A389274983	Apr 21 1996	PT
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 20 2007	CT
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 28 2007	MR
▶	<input type="checkbox"/> ENTERIX	foidh	92485743	Oct 19 2010	PT/SC
▶	<input type="checkbox"/> MACOESSIX	UGVdzg3	987645	Sep 19 1999	PT/SC
▶	<input type="checkbox"/> MANIX	h4873g	0	Nov 1 2003	CT
▶	<input type="checkbox"/> PHENIX	P-000002	A389274983	Mar 21 1996	CT
<input type="checkbox"/>	PHENIX	P-000002	A389274983	Apr 21 1996	PT



Show/Hide filter (by default : Hide)



<input type="checkbox"/>	Patient Name	MRN	Accession #	Study Date	Modality
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 20 2007	CT



Full mode / light mode (by default : Full mode)  
Full mode : Subject - MRN - Accession # - Study date - Modalities - Study description  
Light mode : Subject - Study date - Modalities

All columns can be sort (default : sort by Study date descending)  
The list have no pagination. She load the first 50 studies and if the user scroll at the end, the 50 next studies are loaded.  
When a sort or a filter is apply the list load only the first 50 studies.

## A study with comments

<input type="checkbox"/>	Patient Name	MRN	Accession #	Study Date	Modality
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 20 2007	CT

## A favorite study

<input type="checkbox"/>	Patient Name	MRN	Accession #	Study Date	Modality
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 20 2007	CT

## A linked study

<input type="checkbox"/>	Patient Name	MRN	Accession #	Study Date	Modality
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 20 2007	CT

## A favorite study with

<input type="checkbox"/>	Patient Name	MRN	Accession #	Study Date	Modality
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 20 2007	CT

## A selected study

<input type="checkbox"/>	Patient Name	MRN	Accession #	Study Date	Modality
▶	<input checked="" type="checkbox"/> CERBERIX	uzgcud	3423	Jul 20 2007	CT
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 20 2007	CT

☐ Deselect all selected studies or all series for a study

☒ Selected study

☐ Unselected study

## Mouse over study

<input type="checkbox"/>	Patient Name	MRN	Accession #	Study Date	Modality
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 20 2007	CT
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 20 2007	CT

Clic on a study open it in a viewer

- Add to favorites
- Remove from favorites
- Add a comment (for the first time)
- Show comments
- Get a new link with a
- Get the existing link

## Opened study

<input type="checkbox"/>	Patient Name	MRN	Accession #	Study Date	Modality
▶	<input type="checkbox"/> MACOESSIX	UGVdzg3	987645	Sep 19 1999	PT/SC
<input type="checkbox"/>					
<input type="checkbox"/>					
▶	<input type="checkbox"/> MACOESSIX	UGVdzg3	987645	Sep 19 1999	PT/SC

## A selected series

<input checked="" type="checkbox"/>	
-------------------------------------	--

<input type="checkbox"/>	Patient Name	MRN	Accession #	Study Date	Modality
▶	<input type="checkbox"/> MACOESSIX	UGVdzg3	987645	Sep 19 1999	PT/SC
<input checked="" type="checkbox"/>					
<input type="checkbox"/>					
▶	<input type="checkbox"/> MACOESSIX	UGVdzg3	987645	Sep 19 1999	PT/SC

▼ Opened series

▼ Closed series

## Mouse over series

<input type="checkbox"/>	
--------------------------	--

## Multiple selected studies and series

2 studies and 2 series selected



<input type="checkbox"/>	Patient Name	MRN	Accession #	Study Date	Modality
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 20 2007	CT
▶	<input type="checkbox"/> CERBERIX	uzgcud	3423	Jul 28 2007	MR
▶	<input type="checkbox"/> ENTERIX	foidh	92485743	Oct 19 2010	PT/SC
▶	<input checked="" type="checkbox"/> MACOESSIX	UGVdzg3	987645	Sep 19 1999	PT/SC
▶	<input checked="" type="checkbox"/> MANIX	h4873g	0	Nov 1 2003	CT
<input checked="" type="checkbox"/>					
<input type="checkbox"/>					
▶	<input checked="" type="checkbox"/> MANIX	h4873g	0	Nov 1 2003	CT
▶	<input type="checkbox"/> PHENIX	P-000002	A389274983	Mar 21 1996	CT
<input type="checkbox"/>	PHENIX	P-000002	A389274983	Apr 21 1996	PT

Right clic open a context menu

2 studies and 2 series selected Show the number of studies and series selected




Some as in the context menu (icon and text)

The list is the same for the favorites / recents and inside an album

## Comments


Notifications 

(only in album comment)

 Paul


10/11/2017

ipsum dolor sit amet, consectetur adipisicing elit, sed do eiusmod tempor incididunt ut labore et dolore magna aliqua. Ut enim ad minim veniam, quis nostrud exercitation ullamco laboris nisi ut aliquip ex ea commodo consequat. Duis aute irure dolor in reprehenderit in voluptate velit esse cillum dolore eu fugiat nulla pariatur. Excepteur sint occaecat cupidatat non proident, sunt in culpa qui officia deserunt mollit anim id est laborum.

 Mike


1/10/2017

ipsum dolor sit amet, consectetur adipisicing elit, sed do eiusmod tempor incididunt ut labore et dolore magna aliqua. Ut enim ad minim veniam, quis nostrud exercitation ullamco laboris nisi ut aliquip ex ea commodo consequat. Duis aute irure dolor in reprehenderit in voluptate velit esse cillum dolore eu fugiat nulla pariatur. Excepteur sint occaecat cupidatat non proident, sunt in culpa qui officia deserunt mollit anim id est laborum.

 Mike

Event description

10/10/2017

 Mike (Private message)

[Reply](#)

10/10/2017


ipsum dolor sit amet, consectetur adipisicing elit, sed do eiusmod tempor incididunt ut labore et dolore magna aliqua. Ut enim ad minim veniam, quis nostrud exercitation ullamco laboris nisi ut aliquip ex ea commodo consequat. Duis aute irure dolor in reprehenderit in voluptate velit esse cillum dolore eu fugiat nulla pariatur. Excepteur sint occaecat cupidatat non proident, sunt in culpa qui officia deserunt mollit anim id est laborum.

Use @... for a specific serie  
 Use @... for a specific user  
 Use @... for a specific study (only inside album)

All comments (private and public) are sorted by date the last is down and the first up  
When the user open the comments page, he can see the last comment and scroll up to read the first  
A private message is different  
A private message have a reply button

**Mike** Event description 10/10/2017

Event (only for album)  
Used when a new user join the album or a series is added or removed  
If the event is a new series, a clic open the series


 Use @... for a specific serie  
 Use @... for a specific user  
 Use @... for a specific study (only inside album)

### Sending a public comment

I

Use @... for a specific serie  
Use @... for a specific user  
Use @... for a specific study (only inside album)

Reply to Mike

Cancel

### Sending a private comment after clicking on "reply"

Notifications 

Enable / Disable notification (only for album)

Album list



Name	Study #	Modalities	User#	Messages #	Date	Last event
Album 1	15	CT/MR	5	131	11/12/2017	12/04/2018
Album 2	3	CT/MR/PT	9	18	11/12/2015	12/08/2018
Album 3	27	PT/MR	10	8	11/12/2016	12/05/2018
Album 4	43	CT	5	78	11/12/2017	12/10/2018



Show/Hide filter (by default : Hide)



Name	Study #	Modalities	User#	Comments #	Creation date	Last event
Album 1	15	CT/MR	5	131	11/12/2017	12/04/2018



Full mode / light mode (by default : Full mode)  
Full mode : Name - Study # - Modalities - User # - Comments # - Creation date - Last event  
Light mode : Name - Study # - Modalities - Last event

All columns can be sort (default : sort by Study date descinding)  
The list have no pagination. She load the first 50 albums and if the user scroll at the end, the 50 next albums are loaded.  
When a sort or a filter is apply the list load only the first 50 albums.

A favorite album

Name	Study #	Modalities	User#	Comments #	Creation date	Last event
Album 1 ★	15	CT/MR	5	131	11/12/2017	12/04/2018

Mouse over album

Name	Study #	Modalities	User#	Comments #	Creation date	Last event
Album 1 ☆	15	CT/MR	5	131	11/12/2017	12/04/2018
Album 2	3	CT/MR/PT	9	18	11/12/2015	12/08/2018

Clic on an album open the study list



Add to favorites



Remove from favorites

Album page studies

 Album\_Name ☆

Studies

Comments

Settings

Album description

Lorem ipsum dolor sit amet, erant possit philosophia pro no, saepe utinam tractatos an sit, pri no choro populo. Ea viderer labitur mea, eos clita dolore ut, duo contentiones conclusionemque in. Qui nisl justo ei, molestie salutat

dissentiunt has ex, sea et quis suas. Soluta apeirian usu ei, pri clita aetern

N studies and N series selected

Download


Import in inbox


Delete




Study list

☒

Download

Import to Inbox

Delete

Studies

Comments

Settings

Switch between Study list / Comments / Settings about the current album

The study list is the same like in the inbox

The only difrerence is the context menu when the user do a right clic on the selected studies

Possibilities depend on permissions

Add to favorites

Remove from favorites

Album page settings general


- General
- User
- Token

Album name

Album\_Name 

Album description

Album description    Lorem ipsum dolor sit amet, erant possit philosophia pro no, saepe utinam tractatos an sit, pri no choro populo. Ea viderer labitur mea, eos clita dolore ut, duo contentiones conclusionemque in. Qui nisl justo ei, molestie salutatus dissentiunt has ex, sea et quis suas. Soluta apeirian usu ei, pri clita aeterno dissentias ut. Duo an modo dui, ne vidisse meliore qui.Lorem ipsum dolor sit amet, erant possit philosophia pro no, saepe utinam tractatos an sit, pri no choro populo. Ea viderer labitur mea, eos clita dolore ut, duo contentiones conclusionemque in. Qui nisl justo ei, molestie salutatus dissentiunt has ex, sea et quis suas. Soluta apeirian usu ei, pri clita aeterno dissentias ut. Duo an modo dui, ne vidisse meliore qui.



Notification

- ☐ New Study
- ☒ New comment

 Exit album

- If the user is a admin, he can edit the name and the description.
- A non-admin can just read name and description.
- All user can enable or disable notifications.
- All user can exit the album.
- If the last admin exit the album, a warning message informs the administrator that the album will be admin-free.



# Album page settings user

General

User

Token

Add user

User name

Pierre

Paul

Mike (Admin)

Joe (Admin)

Eric

User settings (only for admin)

Add Studies / Series

Download Studies / Series

Get Studies / Series

Remove Studies /Series

Write comments

Mouse over user

User name

Pierre

Paul

Mike (Admin)

Admin

- Admin

Promote admin
- Delete user

Mouse over admin

User name

Pierre

Paul

Mike (Admin)

An admin can't remove another admin.  
An admin can't downgrade a admin to user

Only an admin can add a new user  
A clic on "add user" open a popup for adding one or more users

Add user

username

username

username

+

Add

Cancel

username

x

For removing an user clic on the "x" when mouse over

User settings (only for admin)

Add Studies / Series

Download Studies / Series

Get Studies / Series

Remove Studies /Series

Write comments

Only an admin can edit the user permissions  
A simple user con only read these permissions

- General
- User
- Token

+ New token

Tokens

☐ Show revoked token

Description	Expiration date	Creation date	Revoke date
Twitter token	12/13/2017	12/13/2016	-
Facebook token	12/13/2017	12/13/2016	-
Twitter token	12/13/2017	12/13/2016	-
Berlin conf token	12/13/2017	12/13/2016	-

Only an admin can see the token setting and generate a new token

Mouse over token

Description	Expiration date	Creation date	Revoke date
Twitter token	12/13/2017	12/13/2016	-
Berlin conf token	12/13/2017	12/13/2016	-

 Revoke the token







☐ Show revoked token

# User Settings

- General
- Token
- Provider

+ New token

Tokens ☐ Show revoked token

Description	Scope	Expiration date	Creation date	Revoke date
Twitter token	 album_name	12/13/2017	12/13/2016	-
Facebook token	 study xxx	12/13/2017	12/13/2016	-
Twitter token	 series yyy	12/13/2017	12/13/2016	-
Berlin conf token	 album_name	12/13/2017	12/13/2016	-
Twitter token	 album_name	12/13/2017	12/13/2016	-
Twitter token	 album_name	12/13/2017	12/13/2016	-

Mouse over token

Description	Scope	Expiration date	Creation date	Revoke date
Twitter token	 album_name	12/13/2017	12/13/2016	-
Facebook token	 study xxx	12/13/2017	12/13/2016	- 

-  token about study
-  token about series
-  token about album
-  revoke token

A clic on the scope open the series / study or album

# New token

Description :

Scope : Album / Series / Study

Album : Album name

Acces permission Read only / Write only / Read-Write

Expiration date / /

Start date / /

Create Cancel

First, the user select the scope and the next comboBox change between Album, series or study.  
The next checkBox use autocomplete.  
Start date an expiration date are optional

## Other

The user can use drag and drop for importing series in her Inbox

When a user do a search in her Inbox, the search is inside the main URL (like in Google)

Internationalization (i18n)