
Header










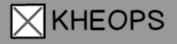

Kheops logo

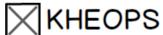

clic: home (Inbox)



Comments notification



Event notification (new study sharing / Job ending / New album)




User (logout / settings)






Notification with number of notification

When the user clic on a notification (event or comment) He will be redirect to the album, comment or studies

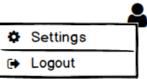

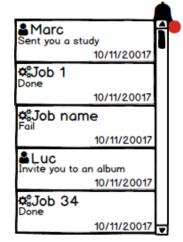



When User clic

When Event notification clic

Sort by date descending Notification source :

User

😋 Job

Notification date Notification

Job

Done

Fail

Album

Paul invite you to the Album_name Paul add studies to the Album_name

Studies

Paul sent you a studies

When Message notification clic

Sort by date descending Notification source :

L User

Album

Notification date

Notification

Message

2018 KHEOPS, Inc.

GitHub Contact About





Kheops logo

GitHub Link


Contact Link


About Link

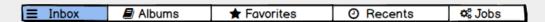

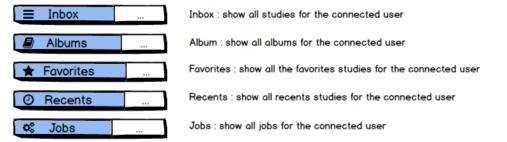

Tabs bar idea from GitHub:

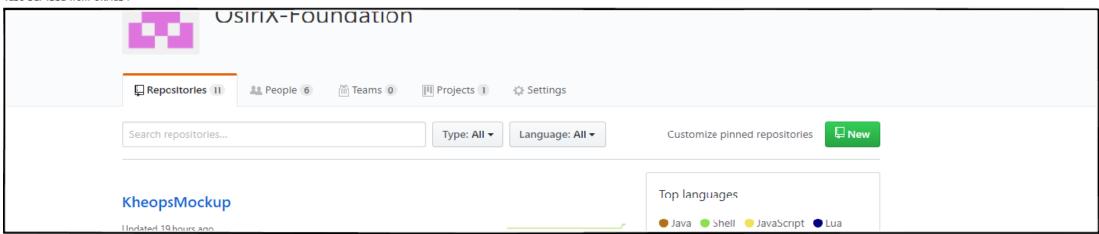

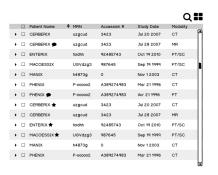




Full mode / light mode (by default : Full mode)
Full mode : Subject - HRN - Accession # - Study date - Modalities - Study descriptio
Light mode : Subject - Study date - Modalities

All columns can be sort (default: sort by Study date descinding)
The list have no pagination. She load the first 50 studies and if the user scroll at the end, the 50 next studies are loaded. When a sort or of liter is apply the list load only the first 50 studies.






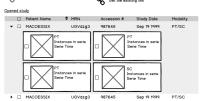



Opened serie
 Closed serie



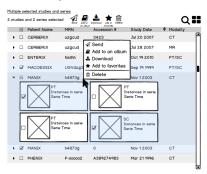



2 studies and 2 series selected Show the number of studies ans series selected



The list is the same for the favorites / recents and inside an album Inside an album, the context menu is defferent



🖴 Paul

10/11/2016

1/10/2017

ipsum dolor sit amet, consectetur adipisicing elit, sed do eiusmod tempor incididunt ut labore et dolore magna aliqua. Ut enim ad minim veniam, quis nostrud exercitation ullamco laboris nisi ut aliquip ex ea commodo consequat. Duis aute irure dolor in reprehenderit in voluptate velit esse cillum dolore eu fugiat nulla pariatur. Excepteur sint occaecat cupidatat non proident, sunt in culpa qui officia deserunt mollit anim id est laborum. Lorem ipsum dolor sit amet. erant possit philosophia pro no. saepe utinam tractatos

≜ Mike

ipsum dolor sit amet, consectetur adipisicing elit, sed do eiusmod tempor incididunt ut labore et dolore magna aliqua. Ut enim ad minim veniam, quis nostrud exercitation ullamco laboris nisi ut aliquip ex ea commodo consequat. Duis aute irure dolor in reprehenderit in voluptate velit esse cillum dolore eu fugiat

Mike

Event description

10/10/2017

Mike (Private message) Reply

10/10/2017

ipsum dolor sit amet, consectetur adipisicing elit, sed do eiusmod tempor incididunt ut labore et dolore magna aliqua. Ut enim ad minim veniam, quis nostrud exercitation ullamco laboris nisi ut aliquip ex ea commodo consequat. Duis aute irure dolor in reprehenderit in voluptate velit esse cillum dolore eu fugiat



Ŧ

Use @... for a specific serie

Use @... for a specific user

Use @... for a specific study (only inside album



All comments (private and public) are sorted by date the last is down and the first up When the user open the comments page, he can see the last comment and scroll up to read the first A private message is different

A private message have a reply button



Event description

10/10/2017

Event (only for album)

Used when a new user join the album or a serie is added or removed

If the event is a new serie, a clic open the serie

Ŧ

Use @... for a specific serie

Use @... for a specific user

Use @... for a specific study (only inside album



Sending a public comment



Use @... for a specific serie

Use @... for a specific user

Use @... for a specific study (only inside album



Sending a private comment after clicking on "reply"






Enable / Disable notification (only for album)



Name	Study #	Modalities	User#	Messages #	Date	Last event
Album 1	15	CT/MR	5	131	11/12/2017	12/04/2018
Album 2	3	CT/MR/PT	9	18	11/12/2015	12/08/2018
Album 3	27	PT/MR	10	8	11/12/2016	12/05/2018
Album 4	43	СТ	5	78	11/12/2017	12/10/2018

Q

Show/Hide filter (by default : Hide)

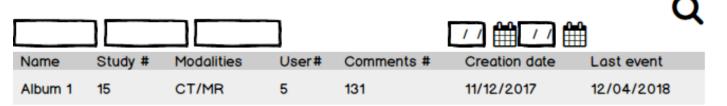



Full mode / light mode (by default : Full mode)

Full mode: Name - Study # - Modalities - User # - Comments # - Creation date - Last event

Light mode: Name - Study # - Modalities - Last event

All columns can be sort (default : sort by Study date descinding)

The list have no pagination. She load the first 50 studies and if the user scroll at the end, the 50 next studies are loaded. When a sort or a filter is apply the list load only the first 50 studies.

A favorite album

Name	Study #	Modalities	User#	Comments #	Creation date	Last event
Album 1	15	CT/MR	5	131	11/12/2017	12/04/2018
Mouse over	album					
Name	Study #	Modalities	User#	Comments #	Creation date	Last event
Name Album 1 ป		Modalities CT/MR	User#	Comments #	Creation date	12/04/2018

Clic on an album open the study list




Add to favorites



Remove from favorites

Album description Lorem ipsum dolor sit amet, erant possit philosophia pro no, saepe utinam tractatos an sit, pri no choro populo. Ea viderer labitur mea, eos clita dolore ut, duo contentiones conclusionemque in. pristi justo ei molestie salutatus dissentiunt has ex sea et quis suas Soluta apeirian usu ei pri clita geterni

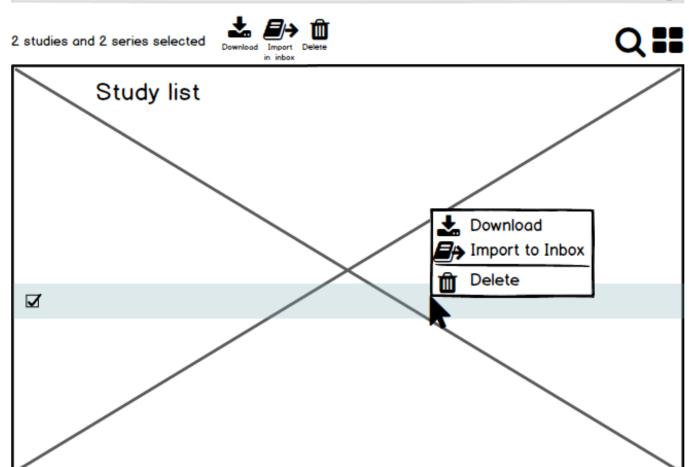

Studies Comments Settings

Switch between Study list / Comments / Settings about the current album

The study list is the same like in the inbox

The only difrerence is the context menu when the use do a right clic on the selected studies

The possibility depand on permissions



